

MethodSCRIPT SDK Example - C



Last revision: March 22 2019

© 2019 PalmSens BV

www.palmsens.com

1.1 Contents:

The example *MethodSCRIPTExample.c* found in the */MethodSCRIPTExample_C* folder demonstrates basic communication with the EmStat Pico using the SDK (C libraries). The example allows the user to start measurements on the EmStat Pico from a Windows PC using a simple C program which makes use of the MethodSCRIPT SDK (C libraries).

1.1.1 Basic Console Example (MethodSCRIPTExample.c)

This example demonstrates how to implement serial communication with the EmStat Pico to

- Establish a connection with the device
- Write a MethodSCRIPT to the device
- Read and parse the measurement data packages from the device

This does not include error handling, method validation etc.

1.2 Communications

The *MSComm.c* from the MethodSCRIPT SDK (C libraries) acts as the communication object to read/write from/to the EmStat Pico.

The functions in the code snippet below are the necessary read/write functions defined in the C example. *WriteToDevice* writes a char to the EmStat Pico and *ReadFromDevice* returns the first byte from the read buffer as soon as there is data available in the buffer.

```
int WriteToDevice(char c)
{
    char writeChar[2] = {c, '\0'};
    if (WriteFile(hCom, writeChar, 1, &dwBytesWritten, NULL))
    {
        return 1;
    }
    return 0;
}

int ReadFromDevice()
{
    char tempChar; //Temporary character used for reading
    DWORD noBytesRead;
    ReadFile(hCom, //Handle of the Serial port
             &tempChar, //Temporary character
             sizeof(tempChar), //Size of TempChar
             &noBytesRead, //Number of bytes read
             NULL);
    return (int)tempChar;
}
```

The read/write functions are required to initiate the *MSComm* communication library.

```
MSComm msComm;
RetCode code = MSCommInit(&msComm, &WriteToDevice, &ReadFromDevice);
```

1.2.1 Connecting to the device

The following code snippet shows how to open a serial com port using the Windows API. In order to use the Windows API for serial communication, the *windows.h* header has to be included in the C program's header file as below.

```
#include <windows.h>
```

A valid handle to the port to which the EmStat Pico is connected is necessary for the read/write functions in the *MSComm* library. The code below can be used to open the com port connected to the device.

```
HANDLE hCom;
hCom = CreateFile(PORT_NAME, GENERIC_READ | GENERIC_WRITE, 0, // must be
opened with exclusive-access
                NULL, // no security attributes
                OPEN_EXISTING, // must use OPEN_EXISTING
                0, // not overlapped I/O
                NULL // hTemplate must be NULL for comm devices
                );

const char* PORT_NAME = "\\.\COM37";
```

The name of the com port connected to the device can be found in the Device Manager in Control Panel in Windows as shown below.

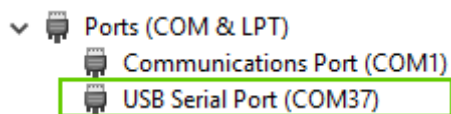


Fig: Available com ports in Device Manager

Once a valid handle is created, the baud rate (230400) for the device has to be set using the device control block object.

1.2.2 Sending the MethodSCRIPT

The methodSCRIPT can be read from a txt file stored in the PC. In this example, the MethodSCRIPT files are stored in the ScriptFiles directory. The code snippet below is used in the example to read the MethodSCRIPT from the file and in turn write it to the device.

```
int SendScriptFile(char* fileName)
{
    FILE *fp;
    char str[100];

    fp = fopen(fileName, "r");
    if (fp == NULL) {
        printf("Could not open file %s", fileName);
        return 1;
    }
    while (fgets(str, 100, fp) != NULL) // Reads a single line from
the script file and writes it on the device.
    {
        WriteStr(&msComm, str);
    }
    fclose(fp);
    return 1;
}
```

1.2.3 Receiving measurement data packages

This example uses the *MSComm* library to receive and parse the data packages from a measurement. In order to read and parse the measurement data packages from the device, the Receive Package function from the *MSComm* library can be used. This function requires a reference to an initiated *MSComm* struct (*msComm*) and it returns the parsed data in the referenced MeasureData struct (*data*)

```
code = ReceivePackage(&msComm, &data);
```

1.2.4 Parsing the measurement data packages

Each measurement data package returned by the function *ReadBuf()* in *MSComm* library, can be parsed further to obtain the actual data values. Here's a set of data packages received from a Linear Sweep Voltammetry (LSV) measurement on a dummy cell with 10 kohm resistance.

```
eM0000\n
Pda7F85F3Fu;ba48D503Dp,10,288\n
Pda7F9234Bu;ba4E2C324p,10,288\n
Pda806EC24u;baAE16C6Dp,10,288\n
Pda807B031u;baB360495p,10,288\n
*\n
\n
```

While parsing a measurement package, various identifiers are used to identify the type of package. For example, In the above sample,

1. 'e' is the confirmation of the "execute MethodSCRIPT" command.
2. 'M' marks the beginning of a measurement loop.
3. 'P' marks the beginning of a measurement data package.
4. "*/\n" marks the end of a measurement loop.
5. "\n" marks the end of the MethodSCRIPT.

The data values to be received from a measurement can be sent through 'pck' commands in the MethodSCRIPT. Most techniques return the data values Potential (set cell potential in V) and Current (measured current in A). These can be sent with the MethodSCRIPT.

In case of Electrochemical Impedance Spectroscopy (EIS) measurements, the following *variable types* can be sent with the MethodSCRIPT and received as measurement data values.

- Frequency (set frequency in Hz)
- Real part of complex Impedance (measured impedance Ohm)
- Imaginary part of complex Impedance (measured impedance in Ohm)

The following metadata values if present can also be obtained from the data packages.

- CurrentStatus (OK, underload, overload, overload warning)
- CurrentRange (the current range in use at the moment)
- Noise (Noise)

1.2.4.1 Parsing the measurement data packages

Each measurement data package begins with the header 'P' and is terminated by a '\n'. The measurement data package can be split into data value packages based on the delimiter ';'. Each of these data value packages can then be parsed separately to get the actual data values.

The type of data in a data package is identified by its variable type:

- The potential readings are identified by the string *da*
- The current readings are identified by the string *ba*
- The frequency readings are identified by the string *dc*
- The real impedance readings are identified by the string *cc*
- The imaginary impedance readings are identified by the string *cd*

For example, in the sample package seen above, the *variable types* are *da7F85F3Fu* - *da* Potential reading and *ba48D503Dp,10,288* - *ba* current reading.

The following 8 characters hold the data value. The data value for the current reading (8 characters) from the above sample package is *48D503Dp*.

The SI unit prefix from the package can be obtained from the parameter value at position 8

In the above sample package, the unit prefix for current data is 'p' which is 1e-12 A.

After obtaining variable type and the data values from the package, the metadata values can be parsed if present.

1.2.4.2 Parsing the metadata values

The metadata values are separated based on the delimiter ',' and each of the values is further parsed to get the actual value.

The first character of each metadata value `metaData[0]` identifies the type of metadata.

'1' – status

'2' – Current range index

'4' - Noise

The status is 1 character hex bit mask. It is converted to int. The status can be obtained as shown in the code snippet below.

For example, in the above sample, the available metadata values for current data are, 10,288. The first metadata value is 10.

1 – metadata status – 0 indicates OK.

The metadata type current range is 2 characters long hex value. If the first bit high (0x80), it indicates a high speed mode current range.

The code below can be used to get current range bits from the package.

The hex value is then converted to int to get the current range string as shown below.

For example, in the above sample, the second metadata available is 288.

2 – indicates the type – current range

88 – indicates the hex value for current range index – 1mA. The first bit 8 implies that it is high speed mode current range.

Sample output:

Here's a sample measurement data package from a LSV measurement on a dummy cell with 10 kOhm resistance and its corresponding output.

Pda7F85F3Fu;ba4BA99F0p,10,288

Output: E (V) = -4.999E-01
 i (A) = -4.999E-01
 Status : OK
 CR : 1mA (High speed)